



Ocean E-VIP Power Height and Tilt

Model no: OCEANEVIP

Designed with the caregiver in mind, the Ocean E-VIP focuses on the safety of the caregiver while providing a secure, non-invasive hygiene experience for the user. The E-VIP's fully powered functions provide smooth, fluid control, through a full range of tilt and height adjustments. Adjustable foot, head and back rests, removable hygiene components, waterproof electronics and stainless steel construction makes care and cleaning for a range of users simple and easy.

Features


- Power adjustable tilting angle -5° to 35° provide a safe seating and repositioning for users with limited weight bearing control and stability
- Seat height is power adjustable between 20" - 41"
- Water-tight power components and stainless steel frame make it possible to clean the complete chair in an institutional washing machine up to 85°C*
- Battery capacity for 40 complete lifting cycles with full load
- Comes complete with adjustable headrest, padded chest belt, collection pan and soft seat
- Compatible with most Ocean/Ocean VIP accessories*

*Except reclining backrest and anti-tippers

Specifications

	Ocean E-VIP OCEANEVIP
Model	
Seat Depth	18"
Seat Width	19"
Seat Height	20" - 41"
Overall Width	29"
Overall Length	40"
Height Under Seat with Pan	13 ½" - 18"
Overall Height (Without Headrest)	42" - 64"
Width Between Armrests	18"
Distance Between Rear Wheels	25"
Maximum Posterior Tilt	35°
Maximum Anterior Tilt	5°
Weight Capacity (lbs.)	330
Product Weight (lbs.)	86

E-VIP Power Height and Tilt Shower Commode

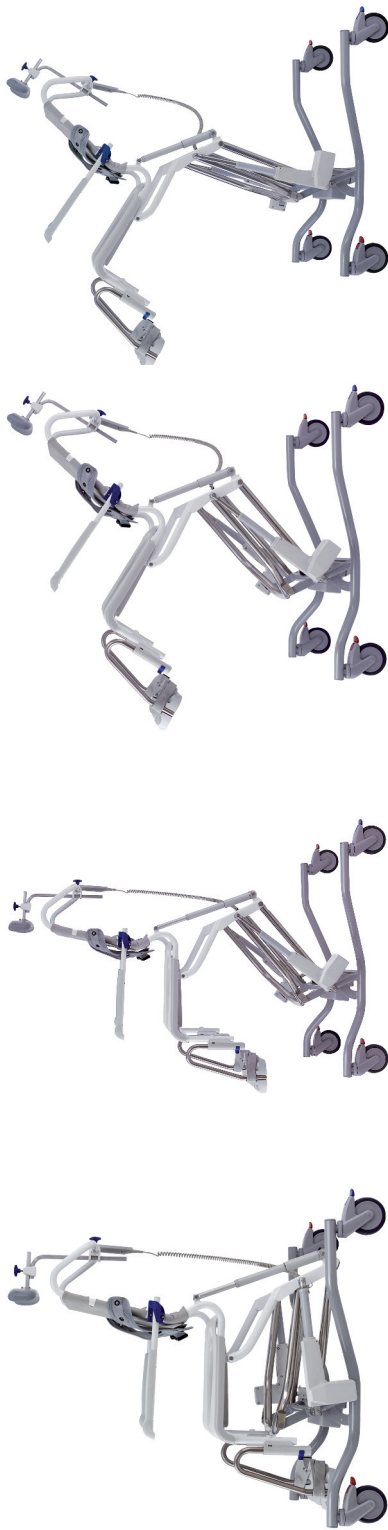

Lifting, transferring and daily hygiene practices are commonly the cause of injury of the caregiver. The Aquatec Ocean E-VIP's height and tilt adjustment enable multiple positioning possibilities. Fully powered operation provides smooth, fluid and precise control, ensuring the safety, comfort and convenience of both the user and caregiver. The E-VIP ensures surface to surface transfer heights are easily accommodated, client re-positioning is smooth and secure, and hygiene and care is performed at safe working height.

Ocean Care Instructions

All components, pan, seat and backrest are safe at temperatures of up to 130°F/85°C for 3 minutes. Aquatec Ocean Series commodes are safe for use with commercial ammonia, chlorinated and alkalis disinfectants and cleaners.

Individual parts can be cleaned using all purpose, commercial cleaning agent detergents using a cloth or brush.

Warranty

Aquatec Ocean shower commodes are covered under a two year limited warranty against manufacturing defects. Manufacturer's warranty does not cover normal wear.

All Aquatec Ocean shower commode accessories are covered for a period of one year.





Yes, you can.®

Invacare Canada LP

www.invacare.ca

570 Matheson Blvd. E.
Unit 8
Mississauga, Ontario
L4Z 4G4 Canada
(800) 668-5324

All rights reserved.
Trademarks are identified
by the symbols TM, SM and [®].
All trademarks are owned
by or licensed to Invacare
Corporation unless
otherwise noted.

Specifications are subject to
change without notification.

© 2018
Invacare Canada
Form No. 08-648C
Rev 09/22