

INVACARE®  
**FULL-ELECTRIC**  
BED



The **Invacare® Full-Electric Bed** offers the greatest convenience for the patient and caregiver. The easy-to use pendant control provides motorized positioning of the upper body and knees, and also adjusts bed frame height for efficient use of accessories and for safe and easy transfer.

Invacare beds, mattresses and rails, excluding therapeutic support surfaces and bariatric beds, when used as a system and sold new after August 1, 2007, meet the voluntary FDA Guidelines for bed rail entrapment. Visit [www.invacare.com/bedsafety](http://www.invacare.com/bedsafety) for additional safety and entrapment information.

## Features

High impact bed end panels are more durable, impact and scratch resistant; washable for easy cleaning

Universal bed ends can be used as either the headboard or footboard on the bed

Rubber caster tread for better tracking with non-marring material

DC low voltage, compact motor system  
Quiet motor operation and fast cycle time on lowering height of bed for efficient bed set-up

Motor system with waterproof enclosures is rated IP66. Washable with standard water pressure. The bed must not be power washed

Color-coded connections on junction box for easier maintenance

Double-insulated junction box provides extra protection against electrical shock

Ergonomic hand pendant comes with durable strain relief and heavy-duty cord. Waterproof casing for easy cleaning, conveniently attaches to bed sheet with clip or rail with hand pendant cradle (sold separately)

Hand pendant cradle available for convenient storage of hand pendant model # 26975

Interchangeability simplifies storage, reduces inventory complexity

## HCPCS Code: E0297



1. Patented universal bed ends
2. Ergonomic hand pendant
3. Color coded connections
4. Emergency manual crank
5. Thermal overload protection
6. Double insulated electrical components/wiring



**Invacare Corporation**  
www.invacare.com

### USA

One Invacare Way  
Elyria, Ohio  
44035-4190  
(800) 333-6900

### Canada

570 Matheson Blvd. E.,  
Unit 8  
Mississauga, Ontario  
L4Z 4G4 Canada  
(800) 668-5324

©2012

Invacare Corporation.  
All rights reserved.  
Trademarks are identified by the symbols ™ and ®. All trademarks are owned by or licensed to Invacare Corporation unless otherwise noted. Specifications are subject to change without notification.  
Form No. 10-123  
Rev. 03/12  
111215



Model #	Dimensions	Sleep Surface	Weight Capacity	Shipping Weight	Carton Dimensions	Regulation Compliance	Limited Warranty
5410IVC	88" L 15"-23" H 36" W	80" L 36" W	450 lb.			UL962 CSA Certified	5 Years Welds 2 Years Mechanical/ Electrical
5000IVC Headspring				45 lb.	36" L 44" W 8" D		
5490IVC Footspring				81 lb.	37" L 47" W 8" D		
5301IVC Universal Bed End				55 lb.	25" L 39" W 7" D		



All references to HCPCS codes: Providers are responsible for determining the appropriate billing codes when submitting claims to the Medicare program, and should consult an attorney or other advisor to discuss specific situations in further detail.

INVACARE®  
**FULL-ELECTRIC  
LOW BED**

---



---

Invacare brings the safety and convenience of a low bed into the Homecare Bed Series with the **Invacare® Full-Electric Low Bed**. The split-spring design incorporated into the standard low bed profile combines the easy, one-person delivery of an Invacare bed with the benefits of a low bed, most often used when injuries from falling out of bed are a concern. The Low Bed ends are fully compatible with existing Invacare full-electric bed frame inventories and the drive shaft, simplifying inventory management for providers.

Invacare beds, mattresses and rails, excluding therapeutic support surfaces and bariatric beds, when used as a system and sold new after August 1, 2007, meet the voluntary FDA Guidelines for bed rail entrapment. Visit [www.invacare.com/bedsafety](http://www.invacare.com/bedsafety) for additional safety and entrapment information.

---

## Features

High impact bed end panels are more durable, impact and scratch resistant; washable for easy cleaning

Universal low bed end can be used as either the headboard or footboard on the Low Bed

Rubber caster tread for better tracking with non-marring material

DC low voltage, compact motor system has quiet motor operation and fast cycle time on lowering height of bed for efficient bed set-up

Washable motor system and frame— Motor system with waterproof enclosures is rated IP66. Washable with standard water pressure. The bed must not be power washed

Color-coded connections on junction box for easier maintenance

Double-insulated junction box provides extra protection against electrical shock

Ergonomic hand pendant comes with durable strain relief and heavy-duty cord. Waterproof casing for easy cleaning, conveniently attaches to bed sheet with clip or bed rail with hand pendant cradle (sold separately)

Hand pendant cradle available for convenient storage of hand pendant (part no. 26975)

Interchangeability simplifies storage, reduces inventory complexity

## HCPCS Code: E0297



1. Bed in lowest position  
2. Bed in highest position



**Invacare Corporation**  
www.invacare.com

**USA**  
One Invacare Way  
Elyria, Ohio  
44035-4190  
(800) 333-6900

**Canada**  
570 Matheson Blvd. E.,  
Unit 8  
Mississauga, Ontario  
L4Z 4G4 Canada  
(800) 668-5324

©2012  
Invacare Corporation.  
All rights reserved.  
Trademarks are identified  
by the symbols ™ and  
®. All trademarks are  
owned by or licensed  
to Invacare Corporation  
unless otherwise noted.  
Specifications are subject  
to change without  
notification.  
Form No. 04-173  
Rev. 03/12  
111218



Model #	Dimensions	Sleep Surface	Weight Capacity	Shipping Weight	Carton Dimensions	Regulation Compliance	Limited Warranty
5410LOW	88" L 9.5"-20" H 36" W	80" L 36" W	450 lb.			UL962 CSA Certified	5 Years Welds 2 Years Mechanical/ Electrical
5000IVC Headspring				45 lb.	36" L 44" W 8" D		
5490IVC Footspring				81 lb.	37" L 47" W 7.5" D		
5301LOW Universal Bed End				55 lb.	25" L 39" W 7" D		

All references to HCPCS codes: Providers are responsible for determining the appropriate billing codes when submitting claims to the Medicare program, and should consult an attorney or other advisor to discuss specific situations in further detail.

INVACARE®  
**SEMI-ELECTRIC**  
BED



The **Invacare® Semi-Electric Bed** combines effortless positioning of the upper body and knees with the economy of manual bed height adjustment. The pendant control provides motorized positioning of the upper body and/or knees. A manual crank located at the foot end of the bed adjusts bed frame height for caregiver assistance.

Invacare beds, mattresses and rails, excluding therapeutic support surfaces and bariatric beds, when used as a system and sold new after August 1, 2007, meet the voluntary FDA Guidelines for bed rail entrapment. Visit [www.invacare.com/bedsafety](http://www.invacare.com/bedsafety) for additional safety and entrapment information.

## Features

High impact bed end panels are more durable, impact and scratch resistant; washable for easy cleaning

Universal bed end can be used as either the headboard or footboard on the bed

Rubber caster tread for better tracking with non-marring material

DC low voltage, compact motor system  
Quiet motor operation and fast cycle time on lowering height of bed for efficient bed set-up

Motor system with waterproof enclosures is rated IP66. Washable with standard water pressure. The bed must not be power washed

Color-coded connections on junction box for easier maintenance

Double-insulated junction box provides extra protection against electrical shock

Ergonomic hand pendant comes with durable strain relief and heavy-duty cord. Waterproof casing for easy cleaning, conveniently attaches to bed sheet with clip or rail with hand pendant cradle (sold separately)

Hand pendant cradle available for convenient storage of hand pendant (part no. 26975)

Interchangeability simplifies storage, reduces inventory complexity

## HCPCS Code: E0295



1. Patented universal bed ends
2. Ergonomic hand pendant
3. Color coded connections
4. Emergency manual crank
5. Thermal overload protection
6. Double insulated electrical components/wiring



**Invacare Corporation**  
www.invacare.com

### USA

One Invacare Way  
Elyria, Ohio  
44035-4190  
(800) 333-6900

### Canada

570 Matheson Blvd. E.,  
Unit 8  
Mississauga, Ontario  
L4Z 4G4 Canada  
(800) 668-5324

©2012

Invacare Corporation.  
All rights reserved.  
Trademarks are identified by the symbols ™ and ®. All trademarks are owned by or licensed to Invacare Corporation unless otherwise noted. Specifications are subject to change without notification.  
Form No. 10-124  
Rev. 03/12  
111215



Model #	Dimensions	Sleep Surface	Weight Capacity	Shipping Weight	Carton Dimensions	Regulation Compliance	Limited Warranty
5310IVC	88" L 15"-23" H 36" W	80" L 36" W	450 lb.			UL962 CSA Certified	5 Years Welds 2 Years Mechanical/ Electrical
5000IVC Headspring				44 lb.	36" L 44" W 8" D		
5890IVC Footspring				77 lb.	37" L 47" W 8" D		
5301IVC Universal Bed End				55 lb.	25" L 39" W 7" D		



All references to HCPCS codes: Providers are responsible for determining the appropriate billing codes when submitting claims to the Medicare program, and should consult an attorney or other advisor to discuss specific situations in further detail.

INVACARE®  
**MANUAL**  
BED



The **Invacare® Manual Bed** offers the same convenient features as electric models, but adjusts manually for greatest economy. Manual cranks located at the foot end of the bed raise and lower upper body and knees, as well as bed height.

Invacare beds, mattresses and rails, excluding therapeutic support surfaces and bariatric beds, when used as a system and sold new and unused after August 1, 2007, meet the voluntary FDA Guidelines for bed rail entrapment. Visit [www.invacare.com/bedsafety](http://www.invacare.com/bedsafety) for additional safety and entrapment information.

## Features

Patented universal bed ends simplify inventory – store just one universal bed end instead of two different bed ends

Patented universal bed ends can be used as either the headboard or footboard on the following Invacare bed models: 5307IVC, 5310IVC, 5410IVC, 5410LOW and VC5310

High-impact bed end panels are stronger than fiberboard and have better impact and scratch resistance

Two locking and two non-locking casters with non-marring rubber tread for secure positioning

Enhanced caster sockets provide better retention during delivery and setup

## HCPCS Code: E0293



1. Patented universal bed ends

**Invacare Corporation**  
www.invacare.com

**USA**  
One Invacare Way  
Elyria, Ohio  
44035-4190  
(800) 333-6900

**Canada**  
570 Matheson Blvd. E.,  
Unit 8  
Mississauga, Ontario  
L4Z 4G4 Canada  
(800) 668-5324

©2015  
Invacare Corporation.  
All rights reserved.  
Trademarks are identified  
by the symbols™ and  
®. All trademarks are  
owned by or licensed  
to Invacare Corporation  
unless otherwise noted.  
Specifications are subject  
to change without  
notification.  
Form No. 11-215  
Rev. 08/15  
1506103



Model #	Dimensions	Sleep Surface	Weight Capacity	Shipping Weight	Carton Dimensions	Regulation Compliance	Limited Warranty
5307IVC	88" L 15"-23" H 36" W	80" L 36" W	Safe Working Load: 450 lb.* Max. Patient Weight: 350 lb.	N/A	N/A	UL962 CSA Certified	5 Years Welds 2 Years Mechanical
5000IVC Headspring				45 lb.	43" L 8" W 36" D		5 Years Welds 2 Years Mechanical
5770IVC Footspring				62 lb.	47" L 9" W 36" D		
5301IVC Universal Bed End				55 lb.	39" L 7" W 24" D		

\*Includes patient, mattress and all accessories.

All references to HCPCS codes: Providers are responsible for determining the appropriate billing codes when submitting claims to the Medicare program, and should consult an attorney or other advisor to discuss specific situations in further detail.



INVACARE®  
**TRAPEZE**  
EQUIPMENT



The **Invacare® Trapeze** unit is an important patient room accessory designed to help the users change positions while in bed and to aid in transferring from bed to chair with easy attendant assistance. Trapeze units have a wide range of height adjustments and handbar positions to maximize patient accessibility.

### Invacare® Trapeze Floor Stand

Model # 7714P

Provides free-standing support for model no. 7740P fixed trapeze unit  
Folds flat for convenient storage and handling

Assisted weight capacity of 250 lb.

### Invacare® Fixed Offset Trapeze Bar with Trapeze

Model # 7740P

Heavy-gauge, octagon steel construction; no tools required for easy assembly

Wall bumpers and padded mounting brackets protect bed ends and wall surfaces

For use on model no. 7714P floor stand or on any Invacare Homecare bed head section

Assisted weight capacity of 250 lb.

### Invacare® Bariatric Floor Stand with Trapeze

Model # BARTRAP

Three locking overhead boom angle positions for comfortable, secure use  
Tool-free design simplifies assembly

Four piece sub-assembly significantly reduces setup time

Lightweight design for easy delivery

Free-standing structure permits use with Bariatric, Home Care or non-medical bed  
Assisted 1000 lb. patient weight capacity

**HCPCS Code: E0940/E0910**



7714P: shown with 7740P

7740P: Includes bar, brackets and grab bar

**HCPCS Code: E0912**



BARTRAP

**Invacare Corporation**  
www.invacare.com

#### USA

One Invacare Way  
Elyria, Ohio  
44035-4190  
(800) 333-6900

#### Canada

570 Matheson Blvd. E.,  
Unit 8  
Mississauga, Ontario  
L4Z 4G4 Canada  
(800) 668-5324

©2012

Invacare Corporation.  
All rights reserved.  
Trademarks are identified by the symbols ™ and ®. All trademarks are owned by or licensed to Invacare Corporation unless otherwise noted. Specifications are subject to change without notification.  
Form No. 11-217  
111223




Model #	Height	Base	Reach	Diameter	Weight	Shipping Weight	Limited Warranty
7714P	38"	40" L 32" W	N/A	1.25"	24 lb.	N/A	3 Years
7740P	61"	N/A	34"	1.25"	17 lb.	N/A	3 Years
BARTRAP	73"	50.5" L 40" W	35.5"	N/A	72 lb.	80 lb.	3 Years


All references to HCPCS codes: Providers are responsible for determining the appropriate billing codes when submitting claims to the Medicare program, and should consult an attorney or other advisor to discuss specific situations in further detail.



# INVACARE® BED ACCESSORIES

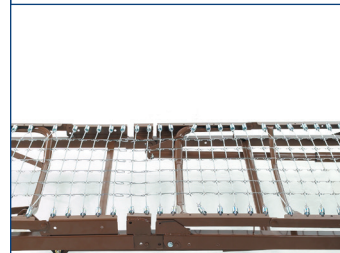

**Invacare® Overbed Tables** are practical, easy-to-use units suitable for both home and institution. The table height easily adjusts for comfort while reading, writing or eating. These tables feature an attractive walnut wood-grained laminate top with a T-molded edge that is designed to resist damage. The heavy-gauge tubular welded steel construction ensures durability and the 1.5" casters ensure smooth mobility of the table in any direction.

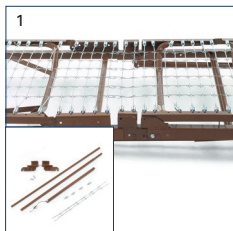

**1. Invacare® Extender Kit**


Model # 5143IVC

Designed to expand 80" bed frame to 84"

Complete with extended head pull tube and drive shaft, link fabric and helicals

Compatible with 5307IVC, 5310IVC, 5410IVC and 5410LOW

Not for use with bariatric beds (BAR600IVC, BAR750), Invacare® Value Care Bed (VC5310) or Invacare® G-Series Bed (G5510)

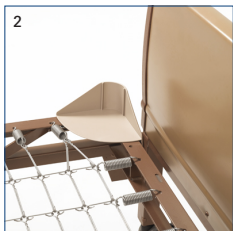

**2. Invacare® Mattress Keepers**

Model # MATTKEEP

Ensures the mattress remains centered on the bed. Maintains mattress position to meet FDA entrapment guidelines. Fits all IVC Series Home Care Beds. Complete kit with plastic mattress keepers and hardware.



**3. Invacare® Tilt-Top Overbed Table**

Model # 6418

Adaptable tilt-top allows three angle positions in either direction using tilt-release lever

Light upward touch

Light downward touch on table top locks it into position



**4. Invacare® Overbed Table**

Model # 6417

Light upward touch on patented spring-loaded locking handle allows table to elevate freely

Light downward touch on table top locks table in position

**Invacare Canada**  
www.invacare.ca

**USA**  
One Invacare Way  
Elyria, Ohio  
44035-4190  
(800) 333-6900

**Canada**  
570 Matheson Blvd. E.,  
Unit 8  
Mississauga, Ontario  
L4Z 4G4 Canada  
(800) 668-5324

©2018  
Invacare Corporation.  
All rights reserved.  
Trademarks are identified by the symbols ™ and ®.  
All trademarks are owned by or licensed to Invacare Corporation unless otherwise noted. Specifications are subject to change without notification.  
Form No. 18-399C  
Rev. 06/18



Model #	Dimensions Table Top	Dimensions Base	Height Adjustment	Diameter	Weight	Weight Capacity	Limited Warranty
6418	30" L 15" W .75" D	26" L 15.5" W	28"-40"	14/16/13 Gauge	22 lb.	25 lb.	1 Year
6417	30" L 15" W .75" D	26" L 15.5" W	28"-40"	14/16/13 Gauge	19 lb.	25 lb.	1 Year

