



INVACARE® MICROAIR® MA600  
THERAPEUTIC SUPPORT SURFACE



Yes, you can.®

**microAIR®**

# Invacare® microAIR® MA600

Alternating Pressure Low Air Loss Mattress System

Medium- High Risk  
Stage I-IV

3-in-1 Alternating Pressure System with  
Low Air Loss aids in the Prevention and  
Treatment of Pressure Injuries



**Highly vapor permeable and oversized pliable quilted nylon top cover** provides low shear, friction and moisture protection



**Cable management system** keeps cords away from the ground to help prevent trip hazards




**Nylon non-slip base cover**  
helps reduce shear and friction during bed articulation






# Features & Benefits

**8 adjustable comfort settings** depending on the patient's therapy requirement

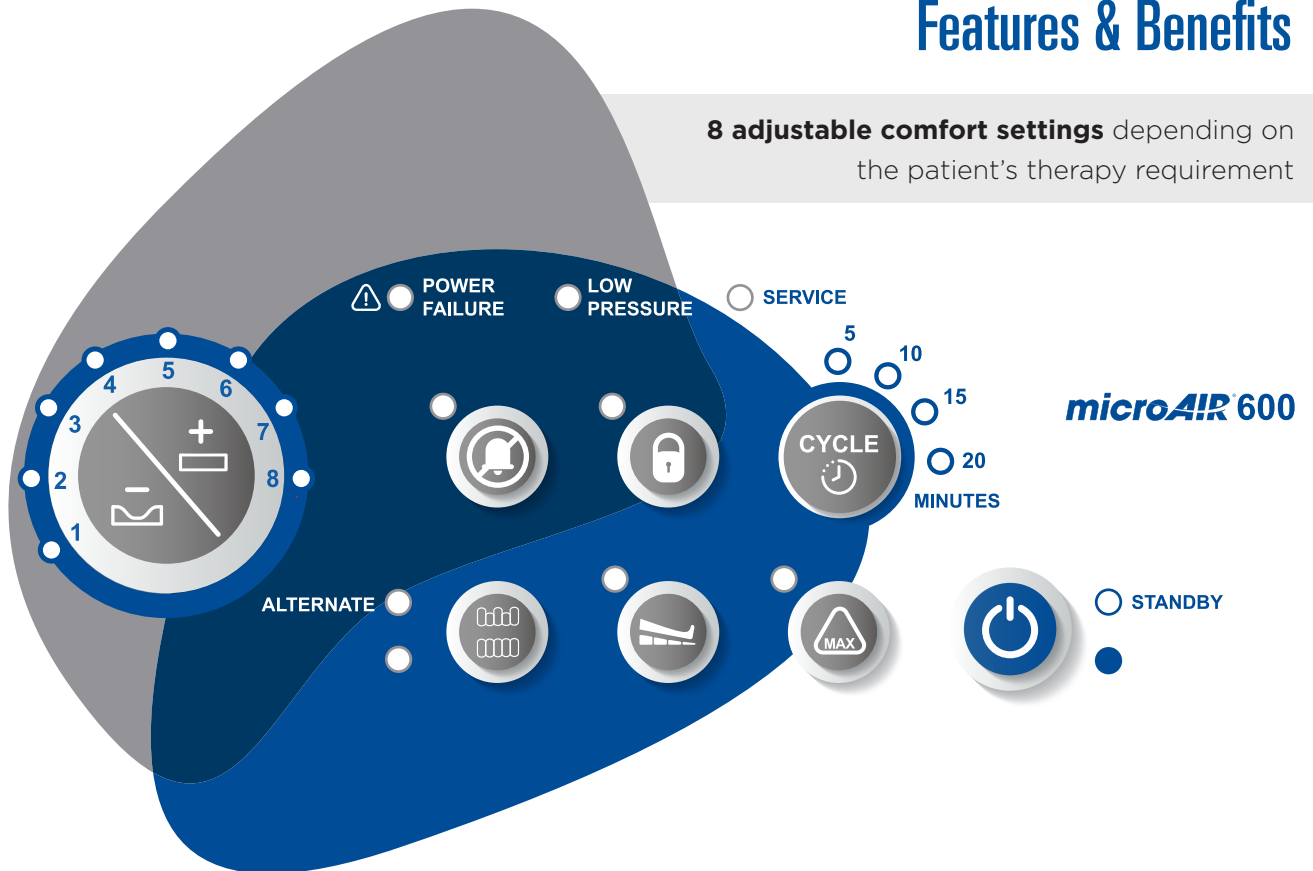

**Advanced 3-in-1 alternating function** inflates and deflates cells in a 3:1 cycle, providing active pressure relief by supporting 2/3 of the body at any given time

**Max Inflate Mode** provides quick maximum inflation within seconds to aid with transfers and nursing procedures

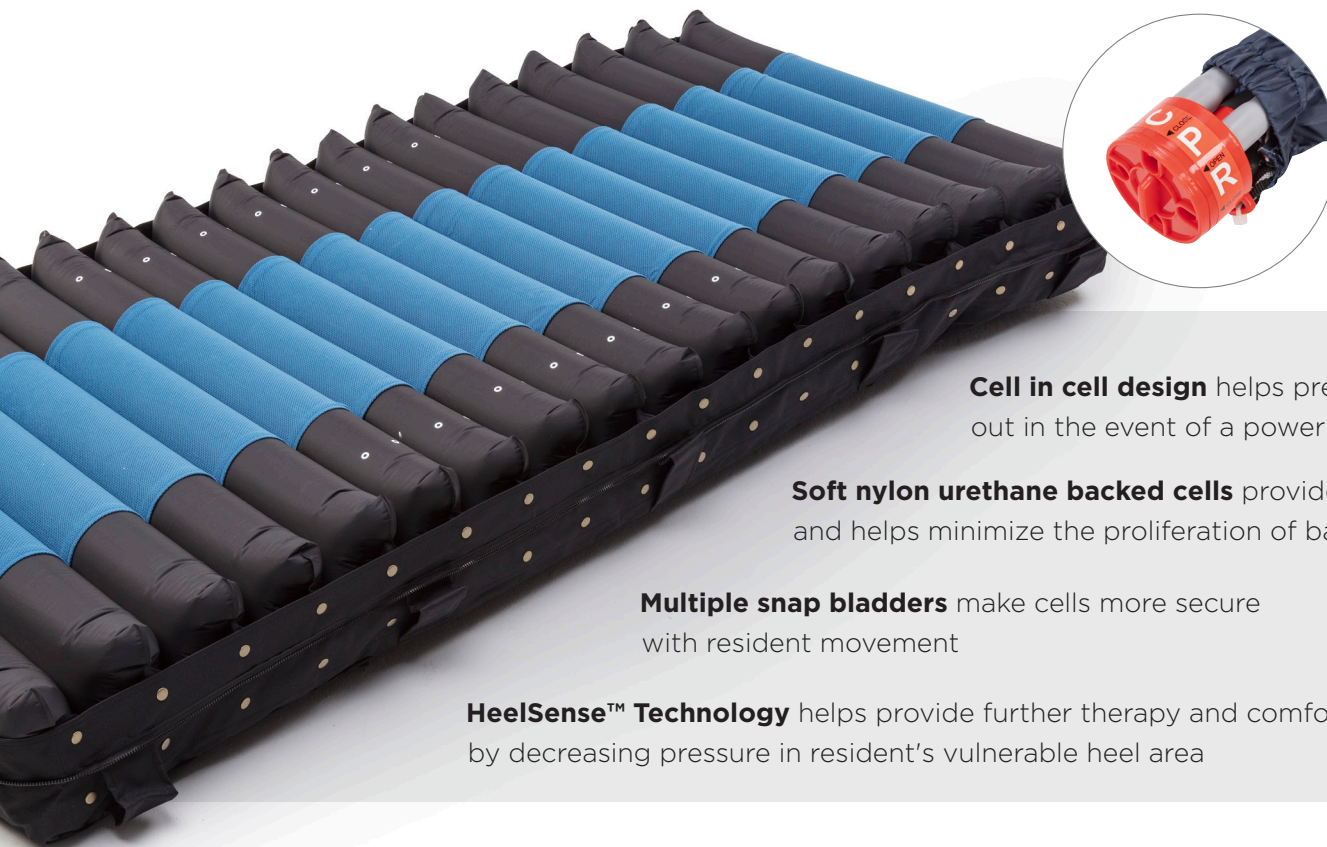

**Cell in cell design** helps prevent bottoming out in the event of a power failure

**Soft nylon urethane backed cells** provide easy cleaning and helps minimize the proliferation of bacterial growth

**Multiple snap bladders** make cells more secure with resident movement

**HeelSense™ Technology** helps provide further therapy and comfort by decreasing pressure in resident's vulnerable heel area

# Request a Quote or Place an Order Today!



**Contact your Inside Sales Representative**

or




**Contact Customer Service at 800.333.6900**

**Invacare Corporation**  
www.invacare.com

**USA**  
One Invacare Way Elyria,  
Ohio  
44035-4190  
(800) 333-6900

**Canada**  
570 Matheson Blvd. E.,  
Unit 8  
Mississauga, Ontario L4Z  
4G4 Canada (800)  
668-5324

## microAIR MA600 Alternating Pressure Low Air Loss Mattress System

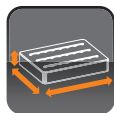

Model #	Dimensions	Weight Capacity	Flow Rate (Pump)	Therapy Type	Alarms	Limited Warranty
MA600	80" L 36" W 8" H	500 lb.	10 LPM	Alternating Pressure Low Air Loss	Audio & Visual	1 Year (Cover/Mattress) 1 Year (Power Unit)
MA600M	microAIR MA600 Mattress only					
MA600P	microAIR MA600 Pump only					

©2019  
Invacare Corporation. All rights reserved. Trademarks are identified by the symbols ™ and ®. All trademarks are owned by or licensed to Invacare Corporation unless otherwise noted. Specifications are subject to change without notification.  
Form No. 19-024  
190203

