



INVACARE® MICROAIR® MA800  
THERAPEUTIC SUPPORT SURFACE



Yes, you can.®

**microAIR®**

# Invacare® microAIR® MA800

Alternating Pressure Low Air Loss Mattress System

**Medium - High Risk  
Stage I-IV**

3-in-1 Alternating Pressure combined with Pulsation and True Low Air is Loss ideal for situations where maceration and moisture management are an issue



**Highly vapor permeable and oversized pliable quilted nylon top cover** provides low shear, friction and moisture protection



**Cable management system** keeps cords away from the ground to help prevent trip hazards



**Nylon non-slip base cover** helps reduce shear and friction during bed articulation







# Features & Benefits

**10 adjustable comfort settings** depending on the patient's therapy requirement

**60 second pulsation** feature helps improve capillary blood flow by continually readjusting to a level surface pressure

**Auto Firm Mode** provides quick maximum inflation within seconds to aid with transfers and nursing procedures

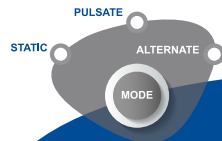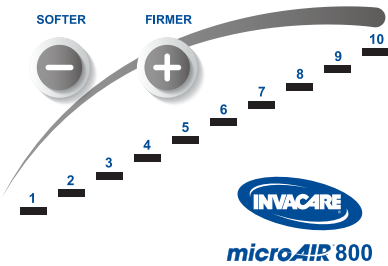

AUTO FIRM LOCK OUT



ALARM RESET



STANDBY

OPERATE

**Advanced 3-in-1 alternating function** inflates and deflates cells in a 3:1 cycle, providing active pressure relief by supporting 2/3 of the body at any given time

**Customizable cycle times** to meet patient's varying needs

**True Low Air Loss Therapy** provides extra comfort for moisture and temperature control



**Foam base** helps prevent bottoming out in the event of a power failure

**Soft nylon urethane backed cells** provide easy cleaning and helps minimize the proliferation of bacterial growth

**External exhaust** keeps the mattress surface cooler and the pump quieter

**HeelSense™ Technology** helps provide further therapy and comfort by decreasing pressure in patient's vulnerable heel area

# Request a Quote or Place an Order Today!



**Contact your Inside Sales Representative**

or




**Contact Customer Service at 800.333.6900**

**Invacare Corporation**  
www.invacare.com

**USA**  
One Invacare Way Elyria,  
Ohio  
44035-4190  
(800) 333-6900

**Canada**  
570 Matheson Blvd. E., Unit  
8  
Mississauga, Ontario L4Z  
4G4 Canada (800) 668-5324

## microAIR MA800 Alternating Pressure Low Air Loss Mattress System

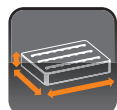

Model #	Dimensions	Weight Capacity	Flow Rate (Pump)	Therapy Type	Alarms	Limited Warranty
MA800	80", 84" L 36", 42", 48" W 10" H	36" 600 lbs. 42", 48" 1000 lbs.	1275 LPM	Alternating Pressure Pulsation True Low Air Loss	Audio & Visual	1 Year (Cover/Mattress) 1 Year (Power Unit)
MA800	microAIR MA800 36" Alternating Pressure Low Air Loss Mattress System					
MA800B42	microAIR MA800 42" Alternating Pressure Low Air Loss Mattress System					
MA800B4284	microAIR MA800 42"W x 84"L Alternating Pressure Low Air Loss Mattress System					
MA800B48	microAIR MA800 48" Alternating Pressure Low Air Loss Mattress System					
MA800M	microAIR MA800 36" Mattress only					
MA800M42	microAIR MA800 42" Mattress only					
MA800M4284	microAIR MA800 42"W x 84" L Mattress only					
MA800M48	microAIR MA800 48" Mattress only					
MA800P	microAIR MA800 Pump only					

©2019 Invacare Corporation. All rights reserved. Trademarks are identified by the symbols ™ and ®. All trademarks are owned by or licensed to Invacare Corporation unless otherwise noted. Specifications are subject to change without notification. Form No. 19-026 190205

